

exposomeShiny User's Guide

Escribà Montagut, Xavier; González, Juan R.

2020-04-02

Contents

1	Overview	5
2	Setup	7
3	Data sets	9
3.1	Exposome dataset	9
3.2	Omics dataset	11
4	Bioconductor packages	13
4.1	rexposome	13
4.2	omicRexposome	13
4.3	CTDquerier	13
5	Analysis flowcharts	15
5.1	Exposome analysis	15
5.2	Exposome-Omic analysis	33
6	General application functionalities	51

Chapter 1

Overview



exposomeShiny is a data analysis toolbox with the following features:

- Data handling: imputation, LOD, transformation, ...
- Exposome characterization
- Exposome-wide association analysis
- Multivariate association
- Omic data integration
- Post-omic data analysis: CTD database

To do so, exposomeShiny relies on previously existent Bioconductor packages (rexposome, omicRexposome and CTDquerier), it uses them in a seamless way so the final user of exposomShiny can perform the same studies that would conduct using the Bioconductor packages but without writing a single line of code.

Chapter 2

Setup

In order to download and setup the environment to launch exposomeShiny, the latest version of the software has to be downloaded from GitHub. To do so open a new RStudio session. Run the following code on the console once the working directory has been setup.


```
# Set working directory
setwd(dir = "/some/path/")

# Download zip
download.file(url = "https://github.com/isglobal-brge/exposomeShiny/archive/master.zip", destfile = "master.zip")

# Unzip the .zip to the working directory
unzip(zipfile = "master.zip")

# Set the working directory inside the downloaded folder
setwd(dir = "/some/path/exposomeShiny-master")
```


Now all the source files are downloaded to the location of chose and the working directory moved to the correct folder, to start the project, open the **Rproj** file by clicking it on the Files explorer of RStudio.

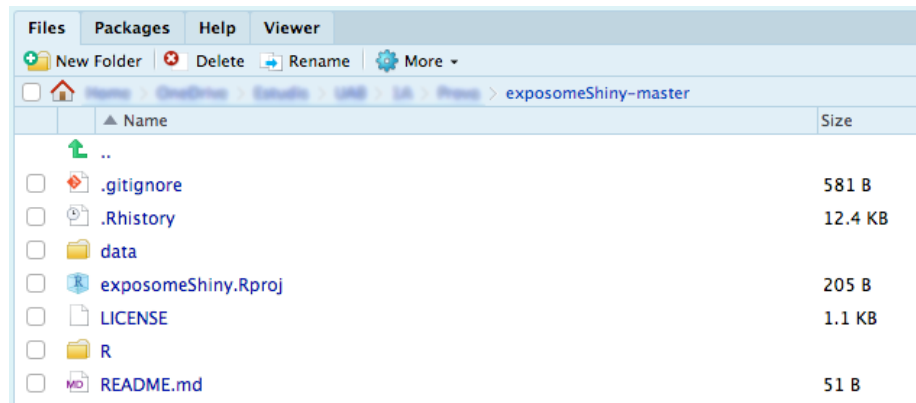

Once the project is loaded, run the following code to install `renv` if it's not present on your R packages library and load the dependencies of `exposomeShiny` on this R session.

```
install.packages("renv")  
  
renv::restore()
```

Now everything is ready to launch the Shiny application. To do so there are two approaches, one is to open the `ui.R` or the `server.R` files that are inside the `R` folder and press `Run App`.

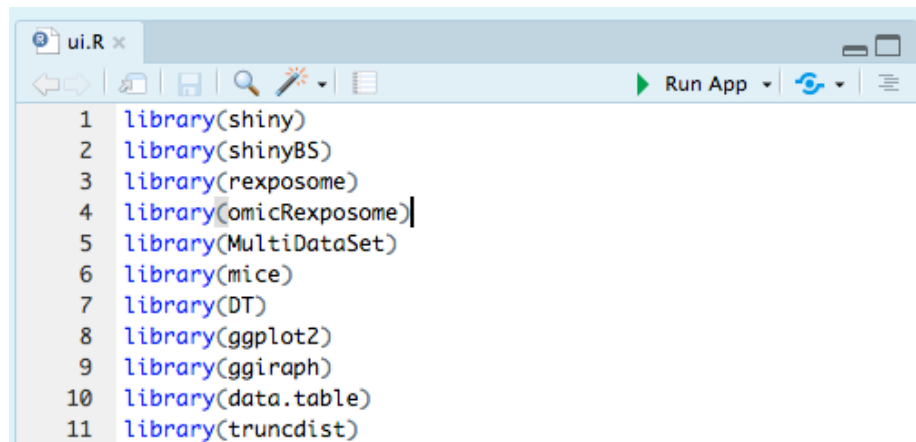

Or the other option is to input the following command on the console.

```
shiny::runApp('R')
```


Chapter 3


Data sets

3.1 Exposome dataset

The exposome is composed of three different files (in `*.csv` format). Those files are referred inside the Shiny as `exposures`, `description` and `phenotypes`. Their content is the following:

- The **exposures** file contains the measures of each exposure for all the individuals included on the analysis. It is a matrix-like file having a row per individual and a column per exposures. It must include a column with the subject's identifier.
- The **description** file contains a row for each exposure and, at last, defined the families of exposures. Usually, this file incorporates a description of the exposures, the matrix where it was obtained and the units of measurement among others.
- The **phenotypes** file contains the covariates to be included in the analysis as well as the health outcomes of interest. It contains a row per individual included in the analysis and a column for each covariate and outcome. Moreover, it must include a column with the individual's identifier.

A visual representation of the three matrices and how they correlate is the following.

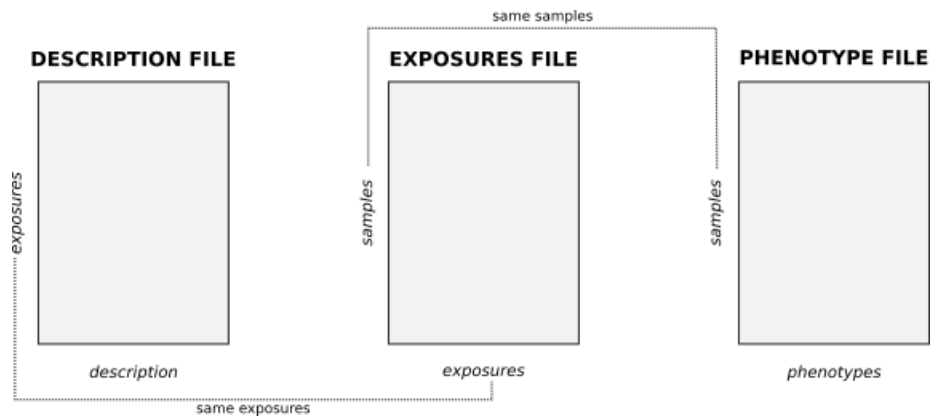

Exposures data file example:

```
id    bde100 bde138 bde209 PFOA    ...
sub01 2.4665 0.7702 1.6866 2.0075 ...
sub02 0.7799 1.4147 1.2907 1.0153 ...
sub03 -1.6583 -0.9851 -0.8902 -0.0806 ...
sub04 -1.0812 -0.6639 -0.2988 -0.4268 ...
sub05 -0.2842 -0.1518 -1.5291 -0.7365 ...
...    ...    ...    ...    ...
```

Description data file example:

exposure	family	matrix	description
bde100	PBDEs	colostrum	BDE 100 - log10
bde138	PBDEs	colostrum	BDE 138 - log10
bde209	PBDEs	colostrum	BDE 209 - log10
PFOA	PFAS	cord blood	PFOA - log10
PFNA	PFAS	cord blood	PFNA - log10
PFOA	PFAS	maternal serum	PFOA - log10
PFNA	PFAS	maternal serum	PFNA - log10
hg	Metals	cord blood	hg - log 10
Co	Metals	urine	Co (creatinine) - log10
Zn	Metals	urine	Zn (creatinine) - log10
Pb	Metals	urine	Pb (creatinine) - log10
THM	Water	---	Average total THM uptake - log10
CHCL3	Water	---	Average Chloroform uptake - log10
BROM	Water	---	Average Brominated THM uptake - log10
NO2	Air	---	NO2 levels whole pregnancy- log10
Ben	Air	---	Benzene levels whole pregnancy- log10

Phenotypes data file example:

```
id    asthma BMI    sex age ...
sub01 control 23.2539 boy 4    ...
sub02 asthma 24.4498 girl 5    ...
```

```
sub03 asthma 15.2356 boy 4 ...
sub04 control 25.1387 girl 4 ...
sub05 control 22.0477 boy 5 ...
...      ...      ...      ...      ...
```

3.2 Omics dataset

The omics data inputted to the Shiny must be provided as an ***.RData**. This file has to contain an `ExpressionSet`, which is an S4 object. This object is a data container of the Bioconductor toolset.

For further information on `ExpressionSet` and how to create and manipulate them, please visit the official documentation and this selected vignette.

Chapter 4

Bioconductor packages

This Shiny application is a front end support for other Bioconductor packages in order to provide a comfortable environment on to conduct different analysis with those packages. In concrete the packages are rexposome, omicRexposome and CTDquerier.

4.1 rexposome

Rexposome is a package that allows to explore the exposome and to perform association analyses between exposures and health outcomes.

4.2 omicRexposome

OmicRexposome is a package that systematizes the association evaluation between exposures and omic data, taking advantage of MultiDataSet for coordinated data management, rexposome for exposome data definition and limma for association testing. Also to perform data integration mixing exposome and omic data using multi co-inherent analysis (omicade4) and multi-canonical correlation analysis (PMA).

4.3 CTDquerier

CTDquerier is a package to retrieve and visualize data from the Comparative Toxicogenomics Database. The downloaded data is formatted as DataFrames for further downstream analyses.


Chapter 5

Analysis flowcharts


5.1 Exposome analysis

As any user would need to do using the Bioconductor packages (rexposome, omicRexposome and CTDquerier) when performing an analysis using an R script, there is some kind of flow (or pipeline) to follow in order to get to the results, this is also true on rexposomeShiny, even though it's a seamless and codeless integration of the packages there's still some need for a flowchart to get the desired results. All the required flowcharts will be detailed with a box flowchart as well as screenshots of exposomeShiny in order to provide extra guidance if needed.


5.1.1 LOD imputation

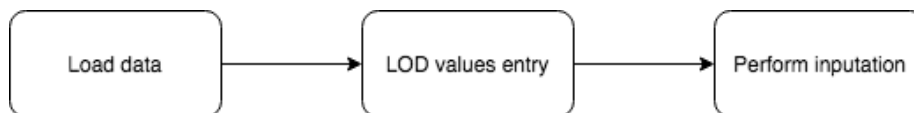


Input the exposures, description and phenotypes files and load them into the application.

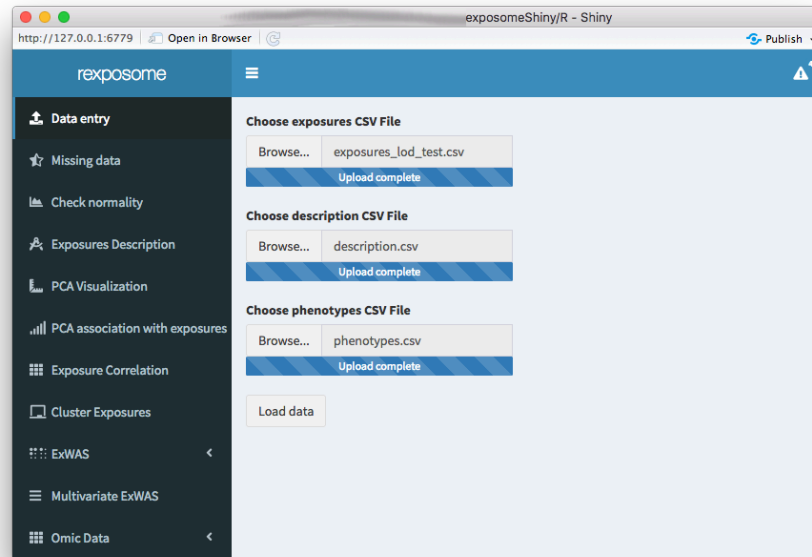


If exposomeShiny detects LODs (limit of detection) on the exposures file (exposures with value: -1), it will prompt the table with the exposures with LOD and double clicking on the desired cell will enable edit mode to input the instrument LOD. There's also the option of selecting "Random imputation" on the imputation method in order to impute with random values instead of $\text{LOD}/\sqrt{2}$.

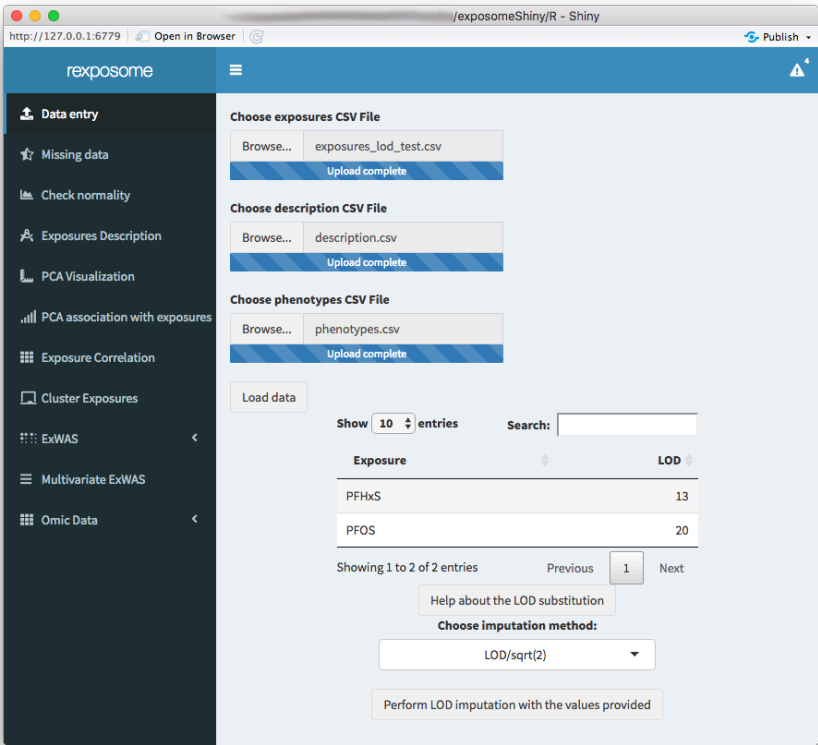

5.1.2 Missing imputation

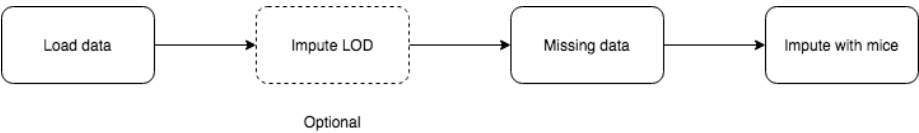


Once the dataset is loaded into the Shiny, look at the “Missing Data” tab to check the percentages of missing data for each exposure present.

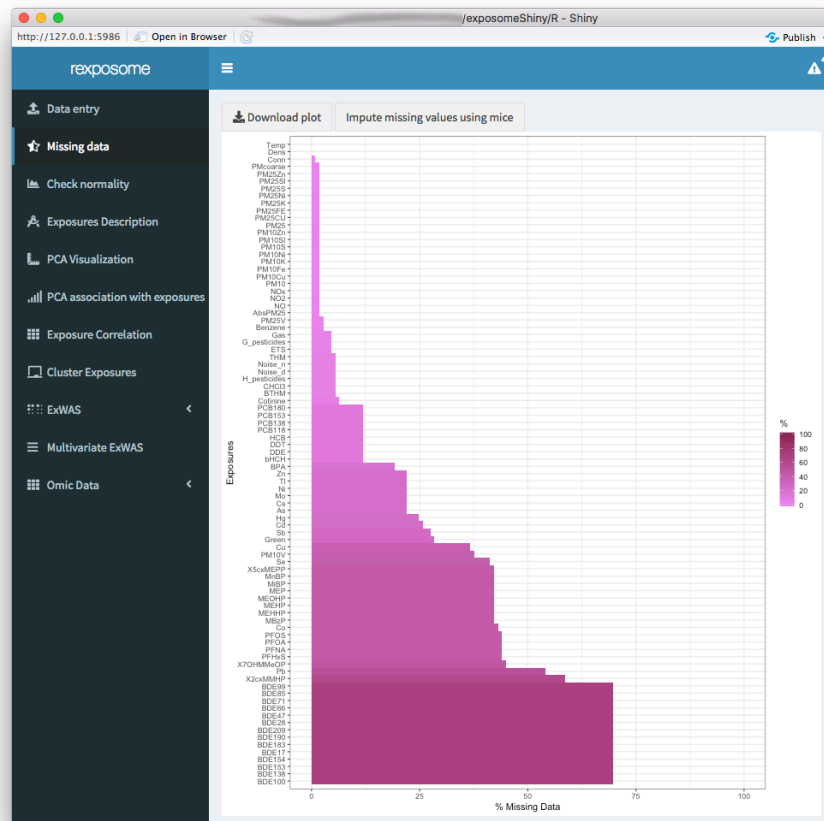


To impute the missing values select “Impute missing values using mice”. After the process finishes, the expect output should be a new missing data graph where there’s no missing of any exposure.

Once the dataset is loaded into the Shiny, look at the “Check Normality” tab to check which exposures are not normal (Normality = false). By selection from the table the desired exposure and clicking the “Plot histogram of selected exposure”, as the label of the button implies, a histogram of the selected exposure from the table can be seen.

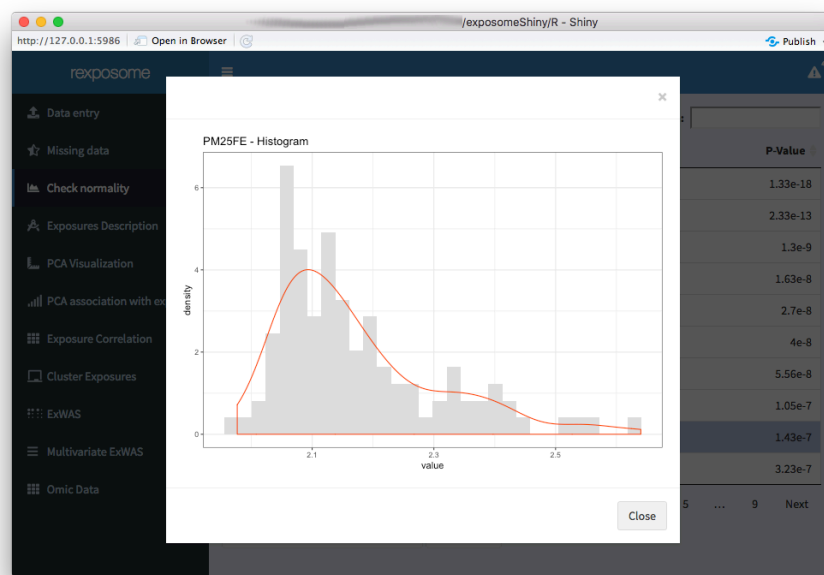


By clicking the “Show false” button, all the non normal exposures are listed with the method that will be applied to normalize, this table can be edited (the “Normalization method” column) by double clicking on the desired row. There are three possible methods to use, “log” (default), “ $\sqrt[3]{}$ ” and “sqrt”. If no method is desired to be applied to an exposure input “none”.

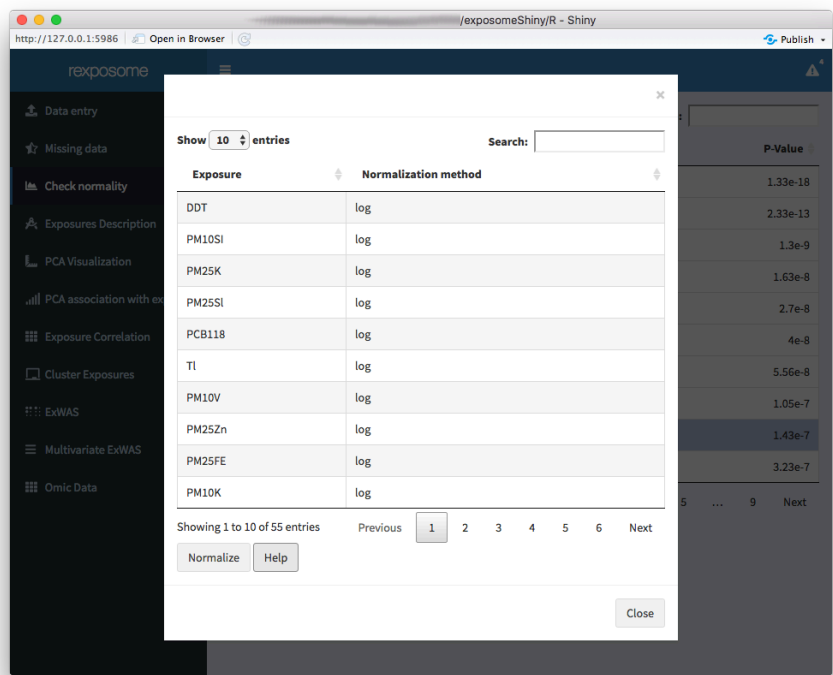


Click “Normalize” and the normalization method selected will be applied, the table on the “Check Normality” tab will be updated with the results of the normalization.

5.1.4 Exposures description

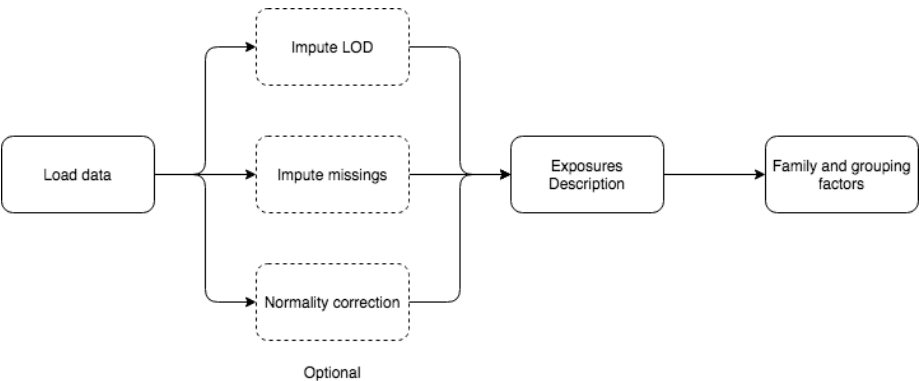


To see all the insights of the exposures dataset loaded into the Shiny, once

loaded it check the exposures description tab, there are three options to dig into the dataset, the family (family of the exposure) to visualize and two grouping factors (phenotypes).

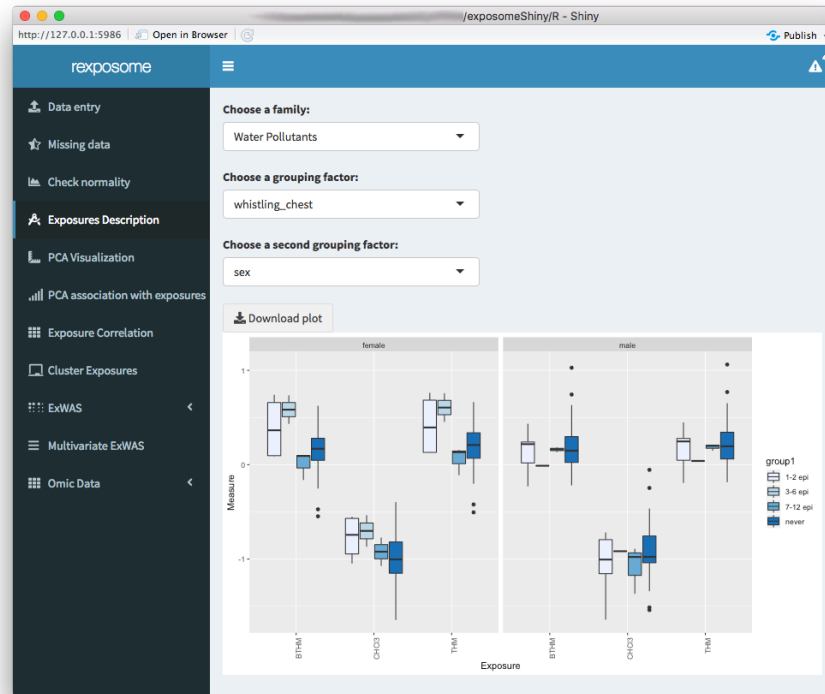


5.1.5 PCA Analysis

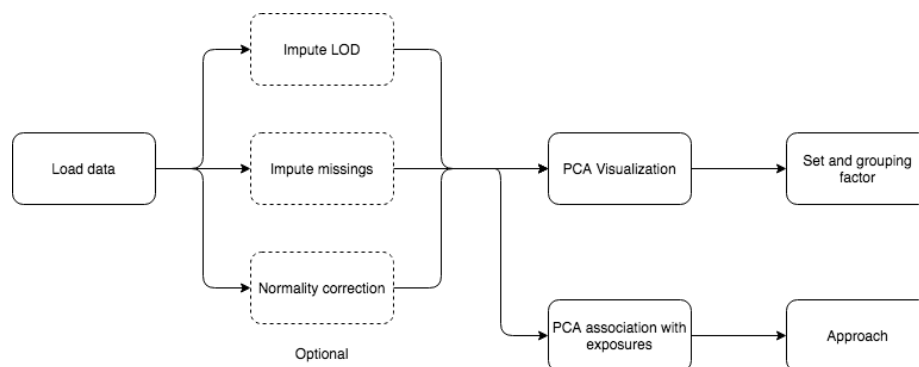


To see the results of a PCA (principal component analysis) study, load the


data and check the PCA Visualization tab, there a set and grouping factor can be choose, it's important noting (as it's already stated on the Shiny) that the grouping parameter only works when the set is selected to “samples”.

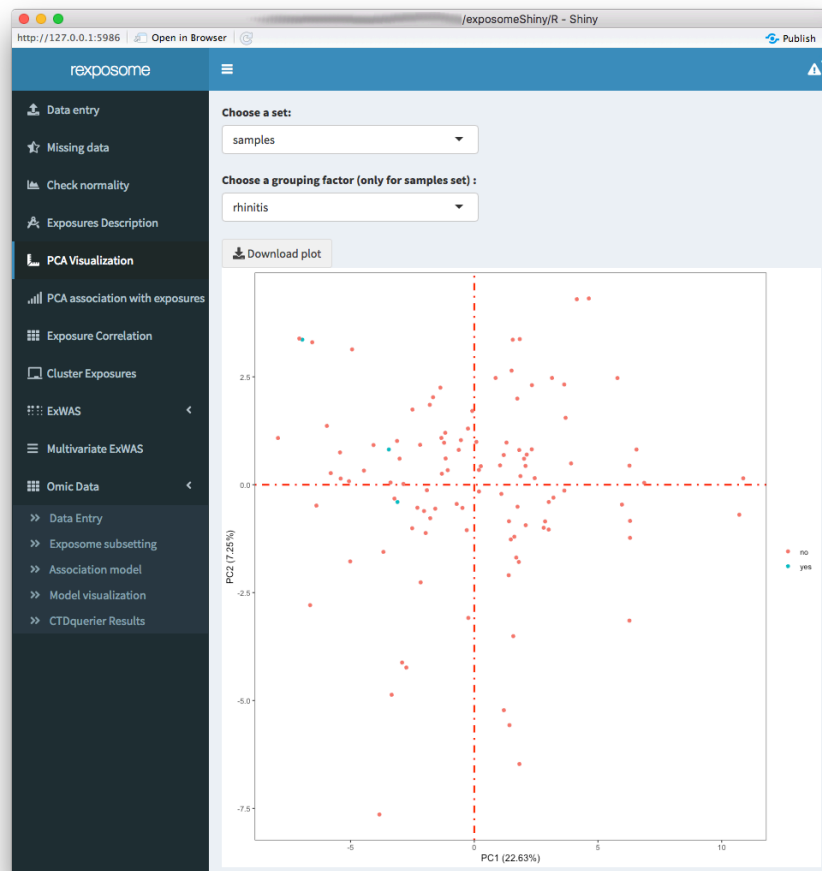

If the association of the PCA analysis with the exposures is desired to visualize, check the “PCA association with exposures” tab, there are two grouping methods to visualize, the phenotypes to principal components and the exposures to principal components.

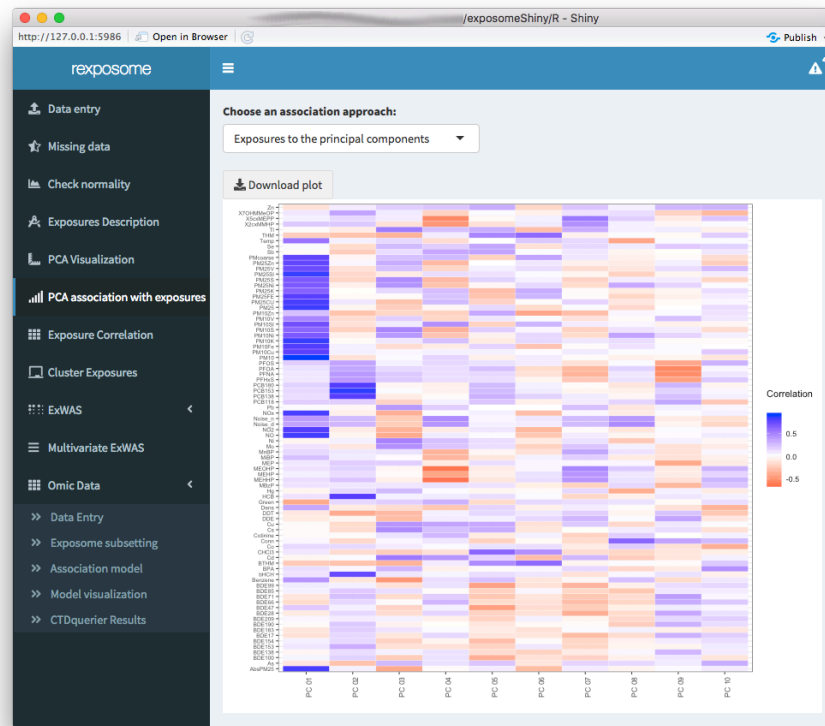


5.1.6 Clusterization and correlation of exposures

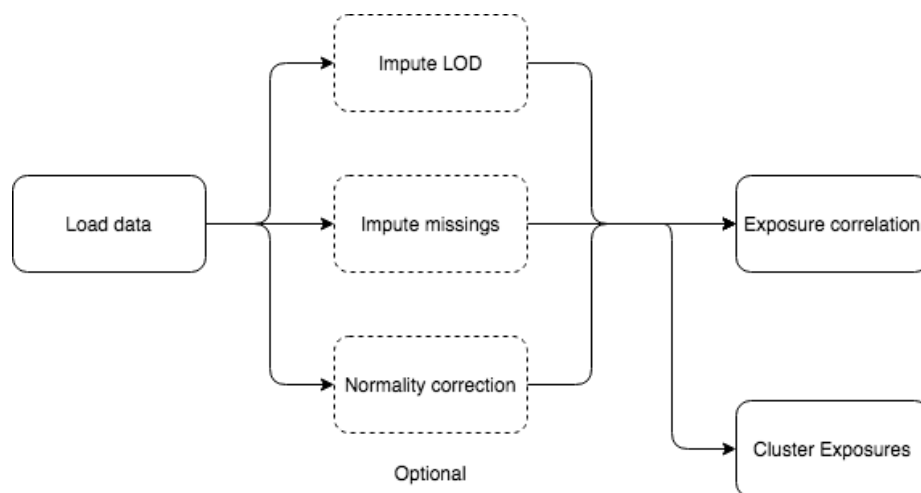


25


[illegible]

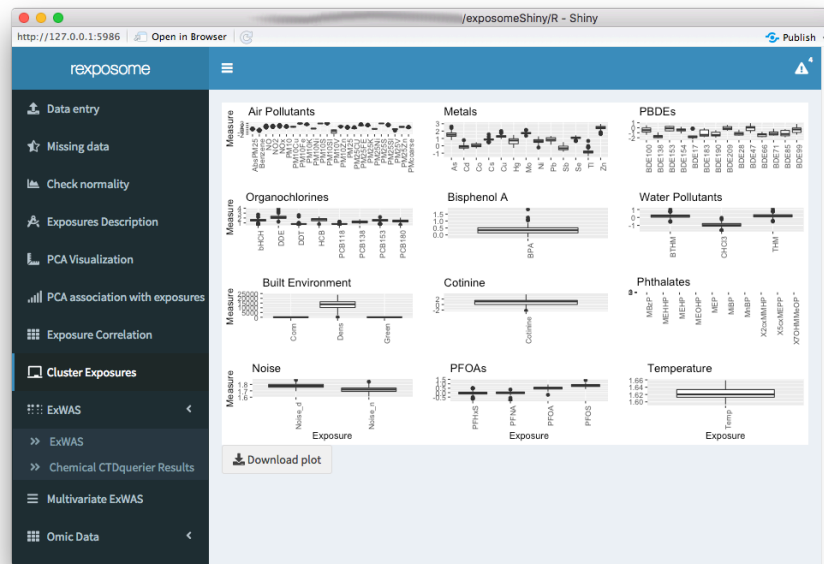

5.1.7 ExWAS

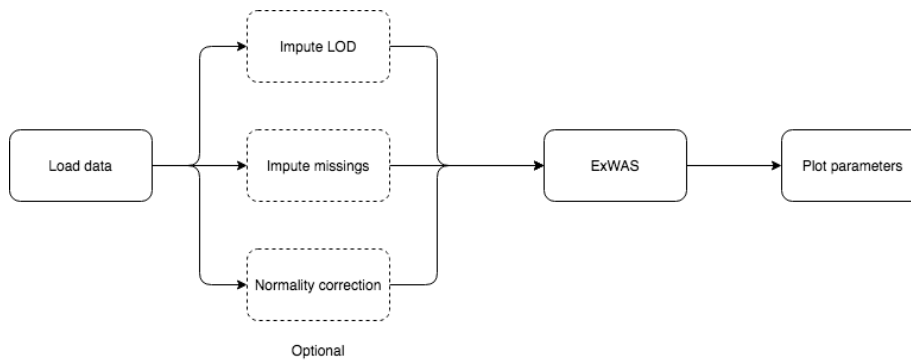


To perform an ExWAS (exposome-wide association) study, check the ExWAS tab and select the adequate parameters for the ExWAS plot, there are two different plot representations, the output variable to choose (phenotype), the output family and as many covariables (phenotypes) as the user wants. There are internal checks to advise the user on which parameters to select depending if the selected outcome is numerical or binomial.

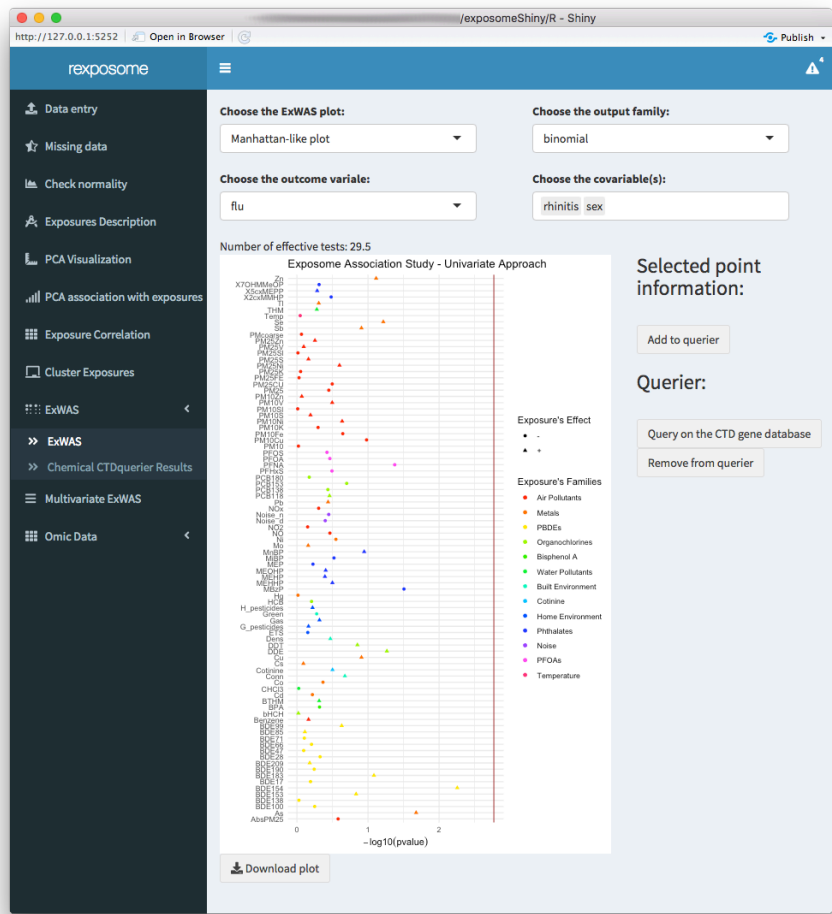



5.1.8 ExWAS - CTDquerier

The ExWAS tab also is able to perform a CTD query of the desired chemicals.

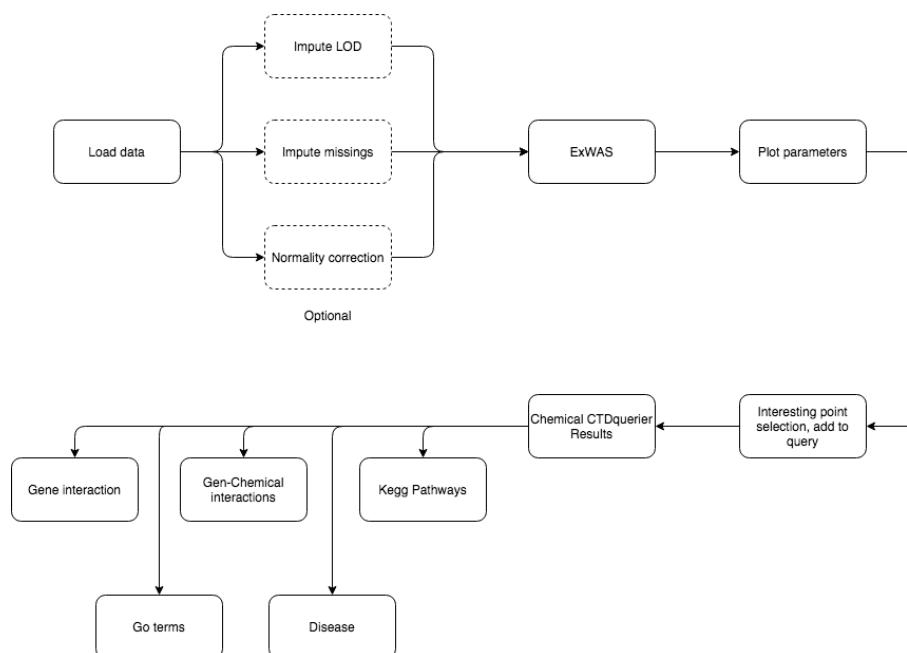


To perform a CTDquerier of chemicals with the results of the ExWAS, click on the desired exposure to preload it into the query, when clicked, a chemical name with it's associated P-Value will appear on the table on the right, if that's the desired chemical to add to the query list click "Add to querier". In the case of adding an unwanted chemical to the query list, select it (or them) by clicking on the Querier list and click on "Remove from querier".

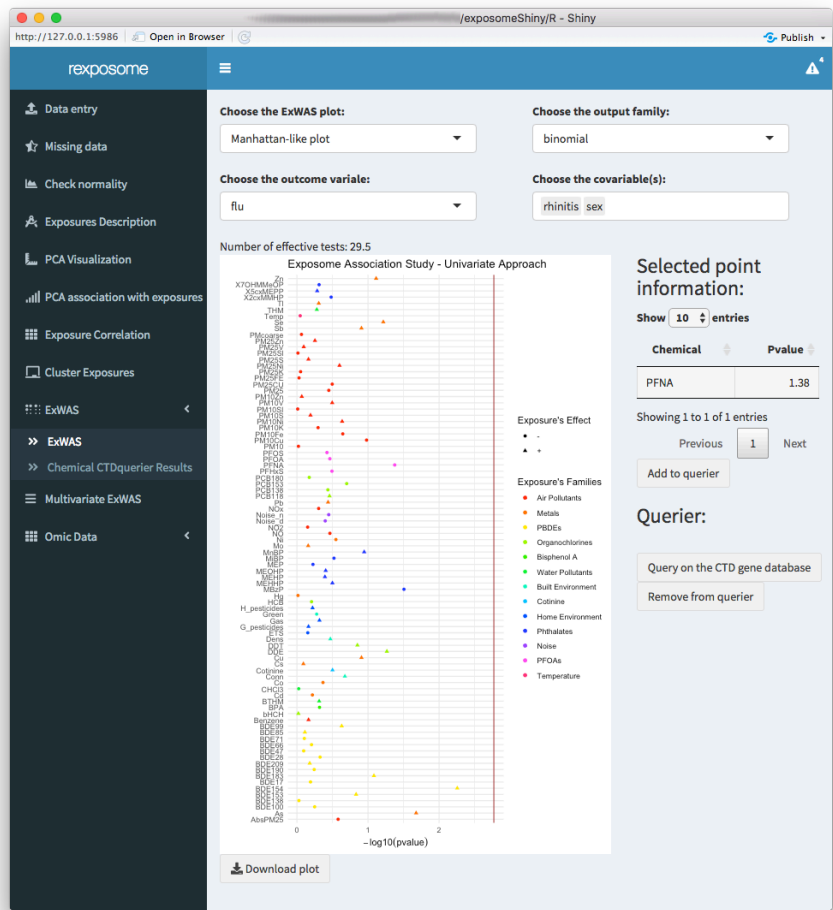


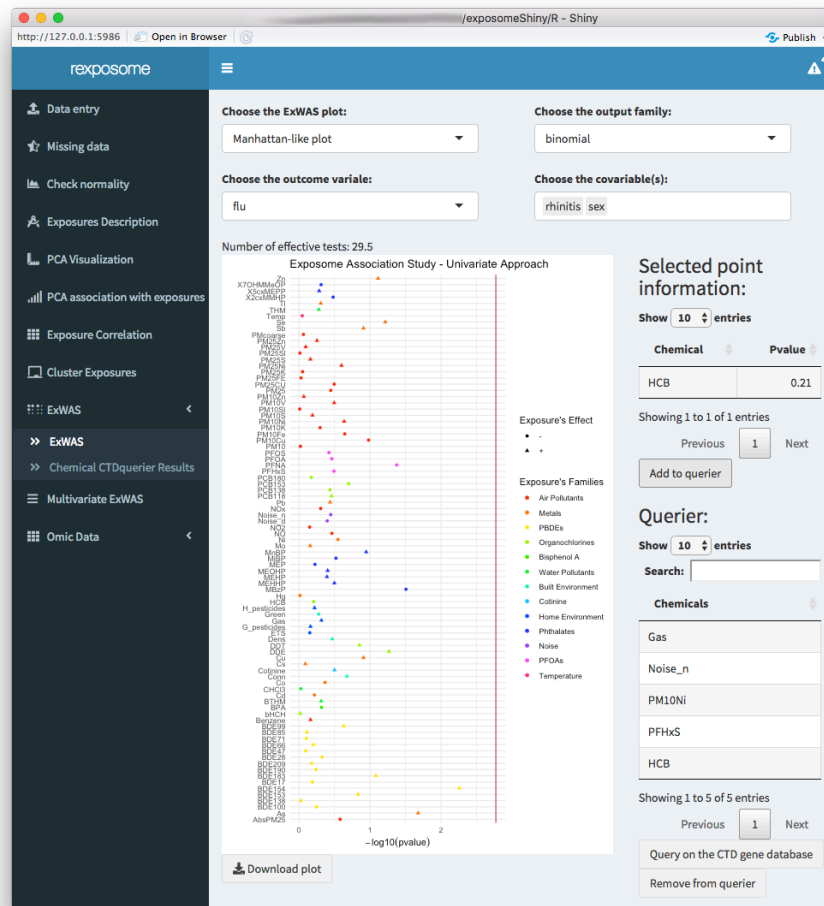


To do the query of the chemicals to de CTD database click on “Query on the CTD gene database” and see the results on the “Chemical CTDquerier Results” subtab. It’s important noting that on the “Kegg pathways” and “Go terms” the input field corresponds to the negative exponent of the filter.

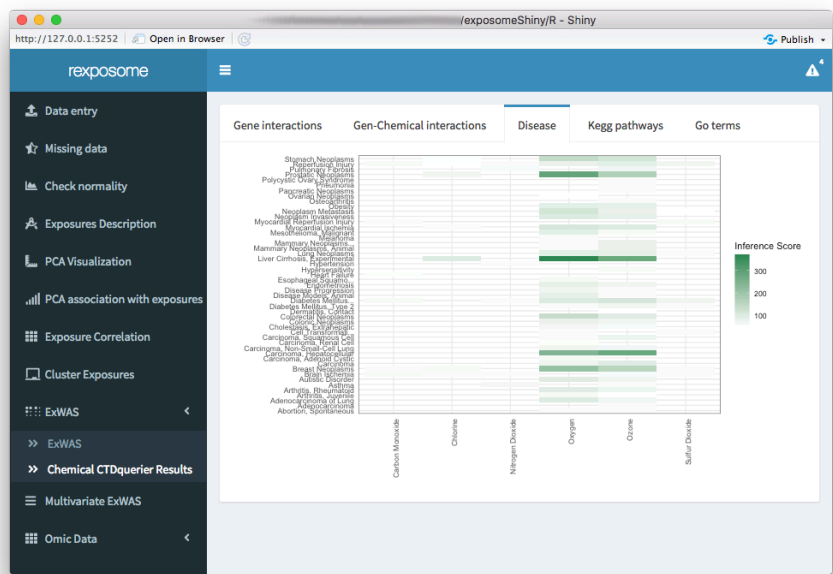


5.1.9 Multivariate ExWAS

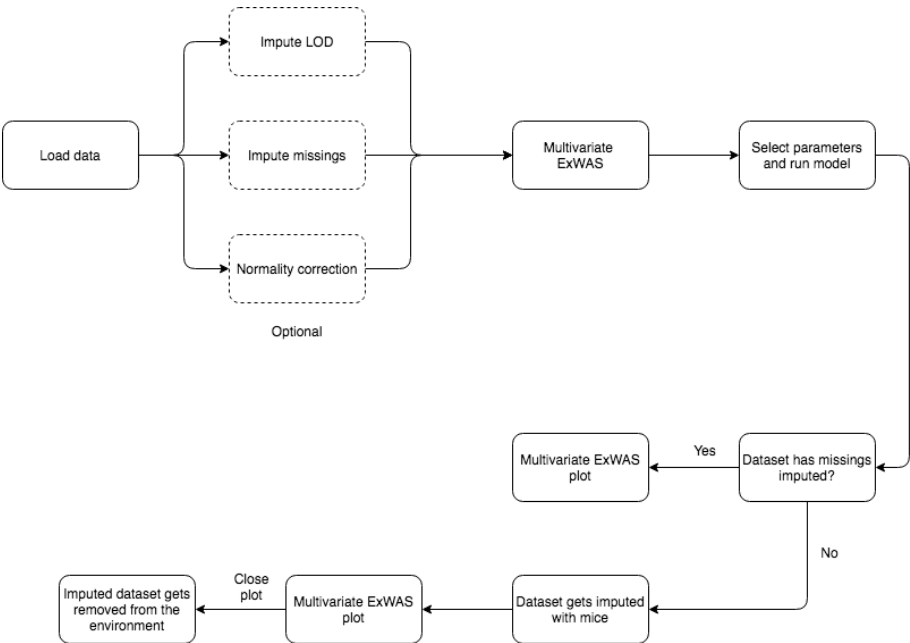

To perform a multivariate ExWAS study, check the Multivariate ExWAS tab and select the desired output parameter, click on run model to generate the plot. As on the ExWAS plot options there's implemented an internal check to advise the user on which parameters to select depending if the selected outcome is numerical or binomial, as the diagrams states if the dataset has not been imputed the missings, it will automatically do it to perform the Multivariate ExWAS, however when closing the plot the imputed dataset will be removed from the environment, so all the other studies performed afterwards will not be altered.

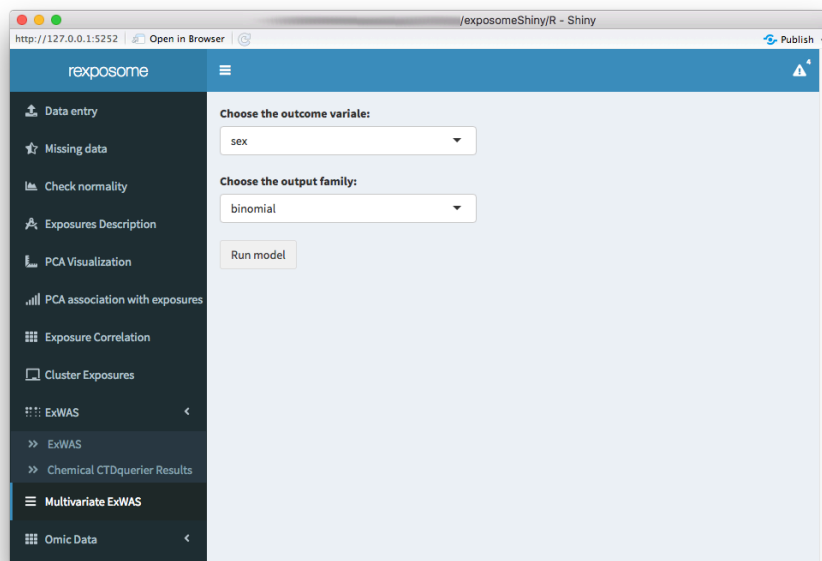

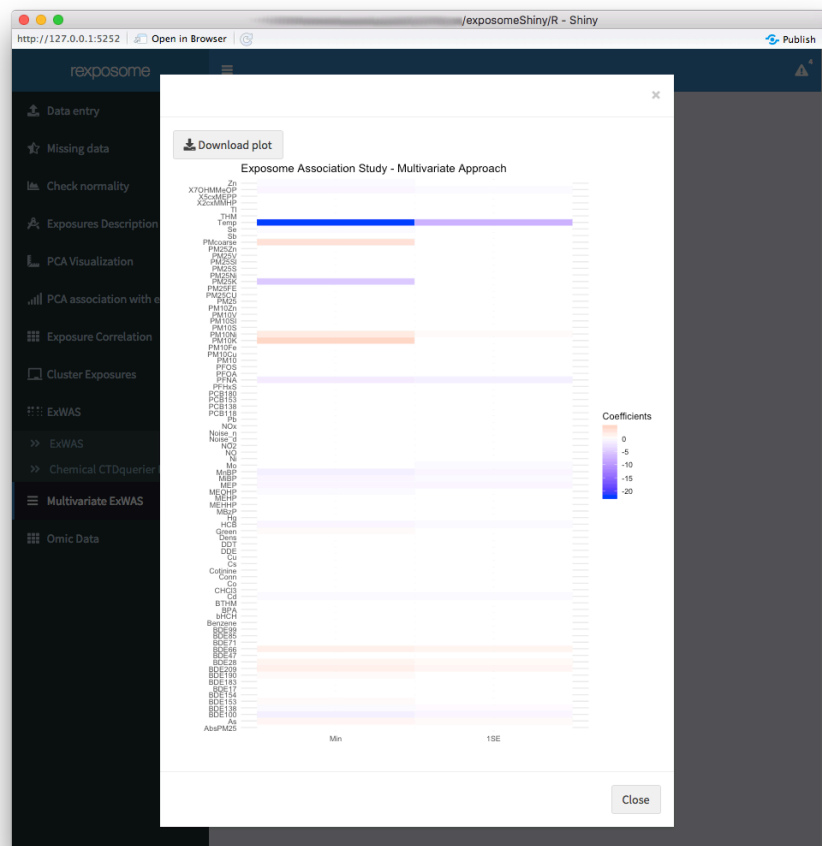


5.2 Exposome-Omic analysis


It's important noting that the maximum size of the omics data is 30 MB, if the omics file to be analyzed is bigger, change the line number 2 of the `server.R` file.

```
# the "30" refers to 30MB, change as needed
options(shiny.maxRequestSize=30*1024^2)
```


5.2.1 Association analysis

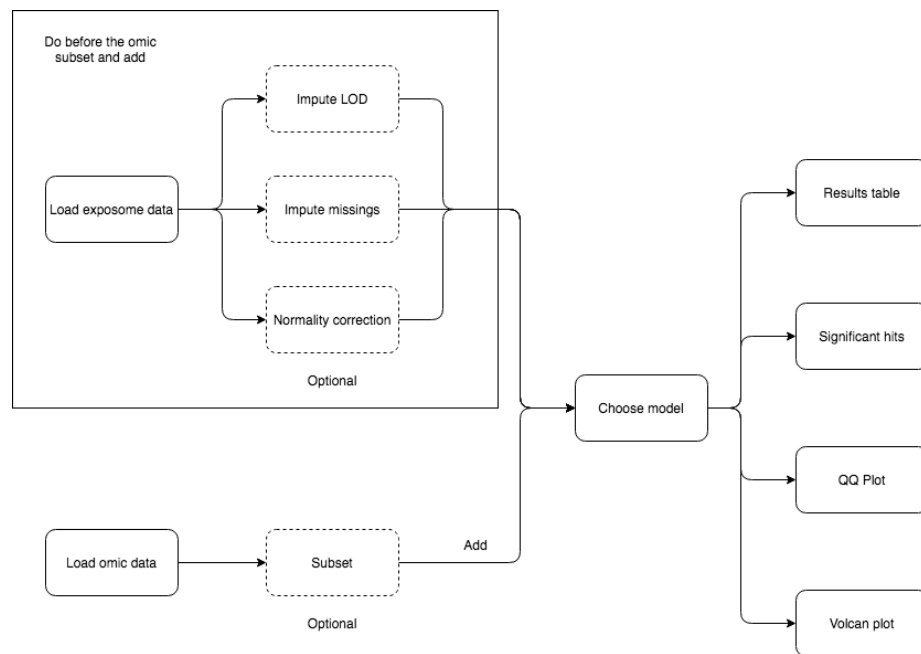


Do first the proceeding of exposome data load and corresponding treatment if desired, then proceed to load the omic dataset on the “Data Entry” subtab of the “Omic Data” tab. The omic data should be provided as a `*.RData` file.

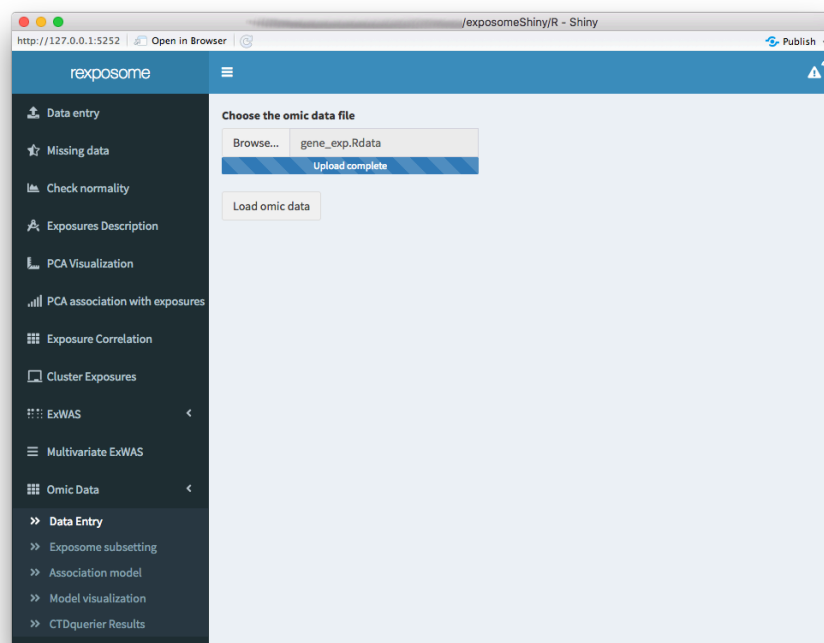

The exposome dataset can be subsetted by families, on the “Exposome subsetting” subtab select the families that are desired to be included in this new set to study, if all the families are desired just don’t input any and proceed to click the “Subset and add”, which will trigger the action to combine the subsetted (or not) exposome dataset with the provided omic dataset.

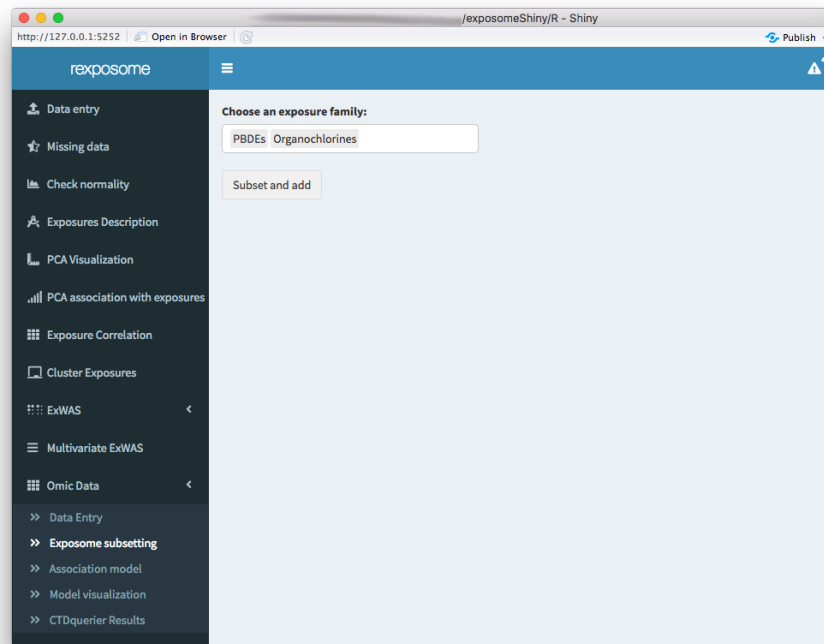

Select the variables for the association analysis and if SVA is wanted on the “Association model” subtab.

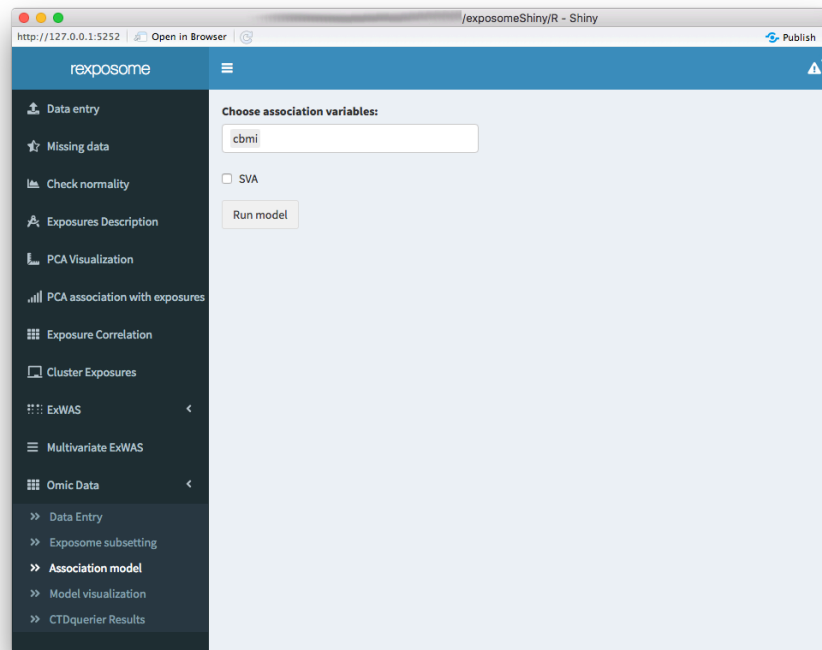


There are tabs to visualize the results of running the association model, all of the are on the “Model visualization” subtab. The “Results table” shows the gene, log of the fold change, p-value and adjusted p-value.


The screenshot shows the reXposome Shiny application interface. The left sidebar contains a navigation menu with options: Data entry, Missing data, Check normality, Exposures Description, PCA Visualization, PCA association with exposures, Exposure Correlation, Cluster Exposures, ExWAS, Multivariate ExWAS, Omic Data, Data Entry, Exposome subsetting, Association model, Model visualization, and CTDquerier Results. The main panel displays the 'Significant hits' tab, which includes a search bar and a table of results. The table has columns for logFC, PValue, and adj.PVal. The results are sorted by logFC in descending order. The table shows 10 entries, with the first entry being TC03002952.hg.1 with a logFC of -1.06, PValue of 0, and adj.PVal of 1. The table also includes a pagination bar at the bottom showing 'Showing 1 to 10 of 20,000 entries' and navigation links for Previous, 1, 2, 3, 4, 5, ..., 2000, and Next.

	logFC	PValue	adj.PVal
TC03002952.hg.1	-1.06	0	1
TC03002591.hg.1	0.84	0	1
TC03000368.hg.1	-1.14	0	1
TC04001638.hg.1	1.07	0	1
TC01004564.hg.1	-1.44	0	1
TC02004830.hg.1	0.86	0	1
TC02002581.hg.1	0.99	0	1
TC04000252.hg.1	0.77	0	1
TC04001957.hg.1	0.68	0	1
TC03001355.hg.1	-2.55	0	1


The “Significant hits” shows the exposure, hits and lambda.

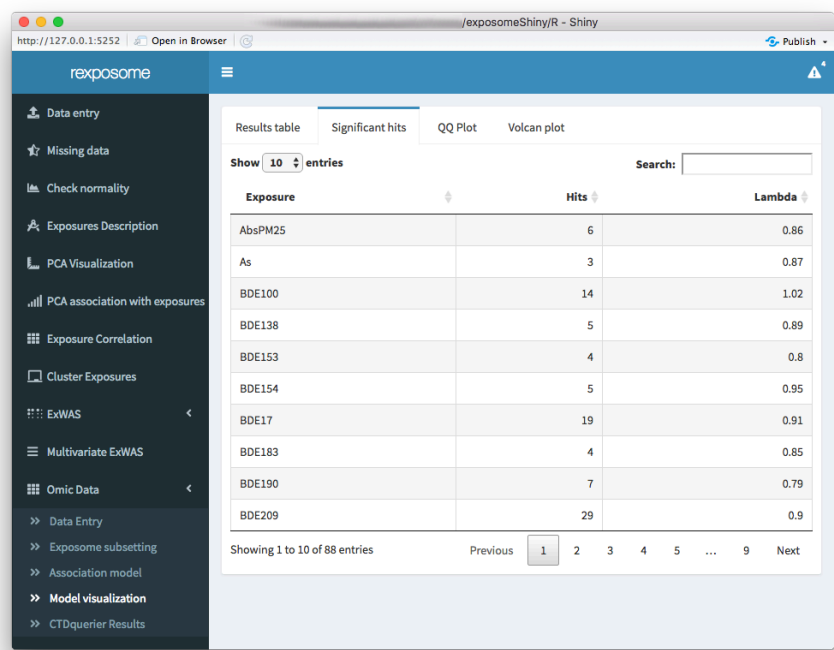


The QQ Plot shows a QQ plot (expected vs. observed $-\log_{10}(\text{p-value})$) for the selected exposure.

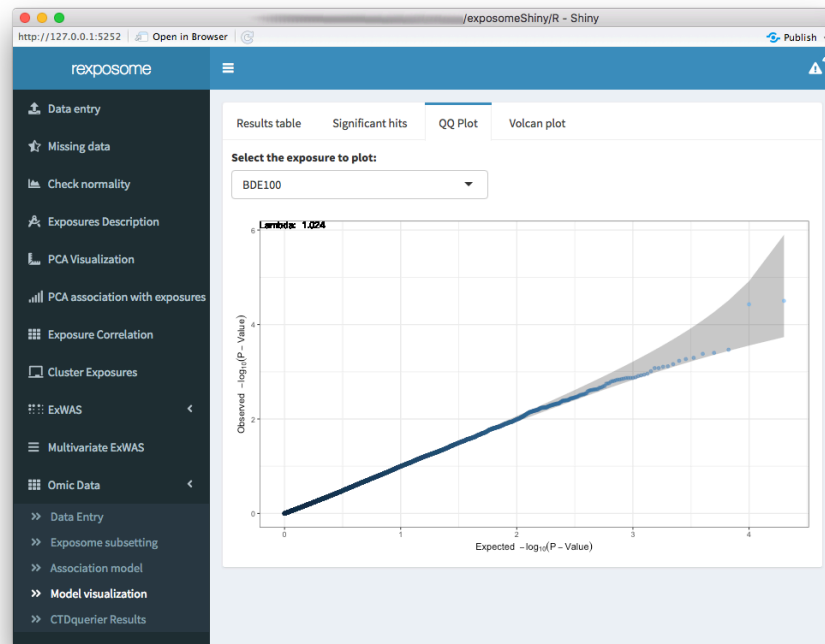


The Volcan plot shows a volcano plot ($\log_2(\text{fold change})$ vs $-\log_{10}(\text{p-value})$). For this plot there are two input cells to adjust the horizontal and vertical limit lines to filter out the results.

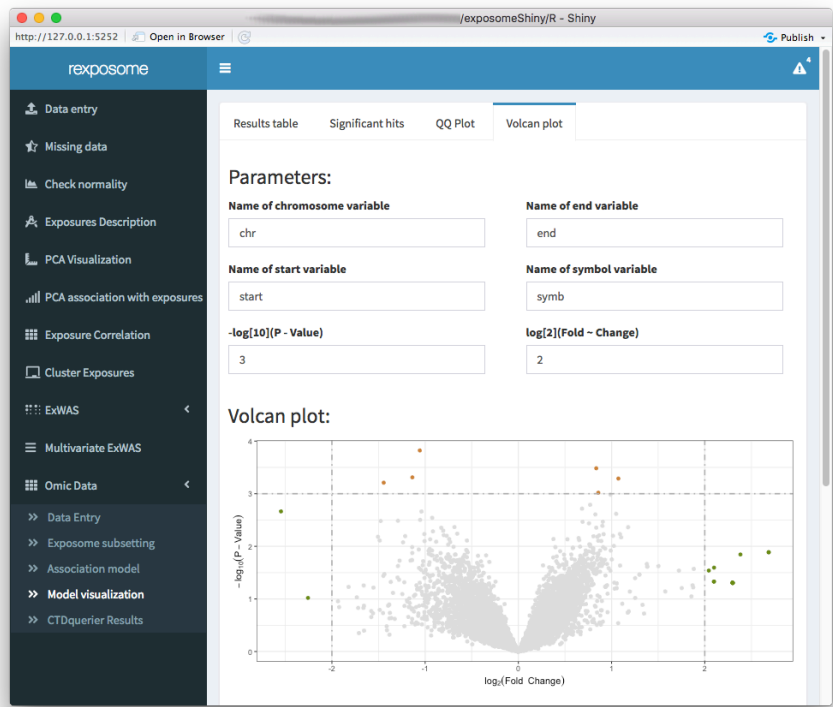


5.2.2 CTD querier

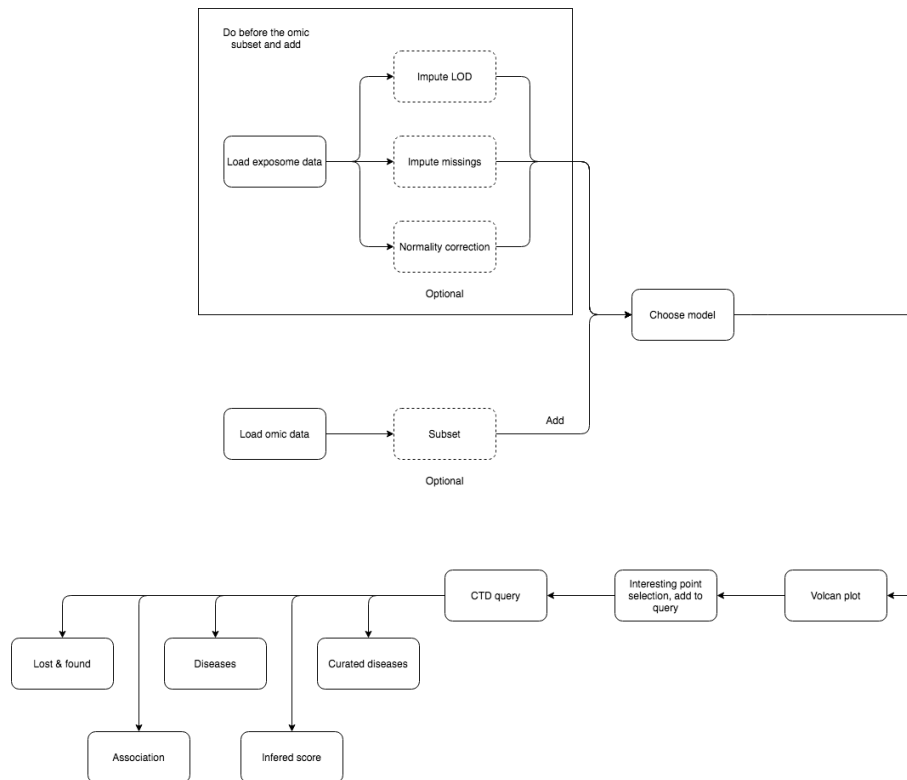


To perform a CTD querier study of the exposome-omic analysis, as before, load both datasets and run the desired model with them, check the Volcan plot and click on the desired point on the Volcan plot, the information of the selected plot will appear on the table below the plot (sometimes there are many points close so more than one rows can appear on the table), select from the table the desired point to add to the query and click “Add to querier”. It’s important noting that when trying to add to the querier the Shiny will find on the fields of the omic dataset that the user specifies on top of the plot. If the search does not return any symbol a prompt will appear, however if it’s found it will be added to the lower table corresponding to the genes to query.

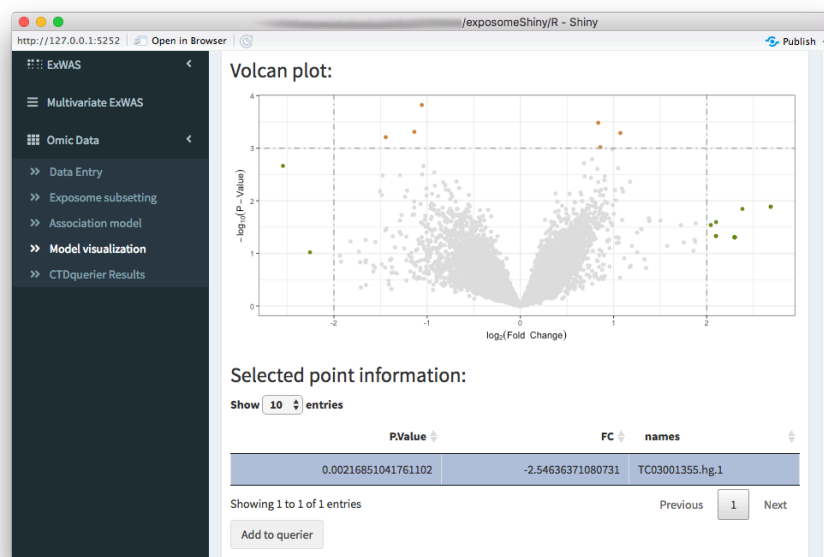

If by mistake some gene (or genes) were introduced to the querier, select them by clicking on the table row and click “Remove from querier”. Click on “Query selected genes on the CTD gene database” to perform the query of all the symbols of the querier list.

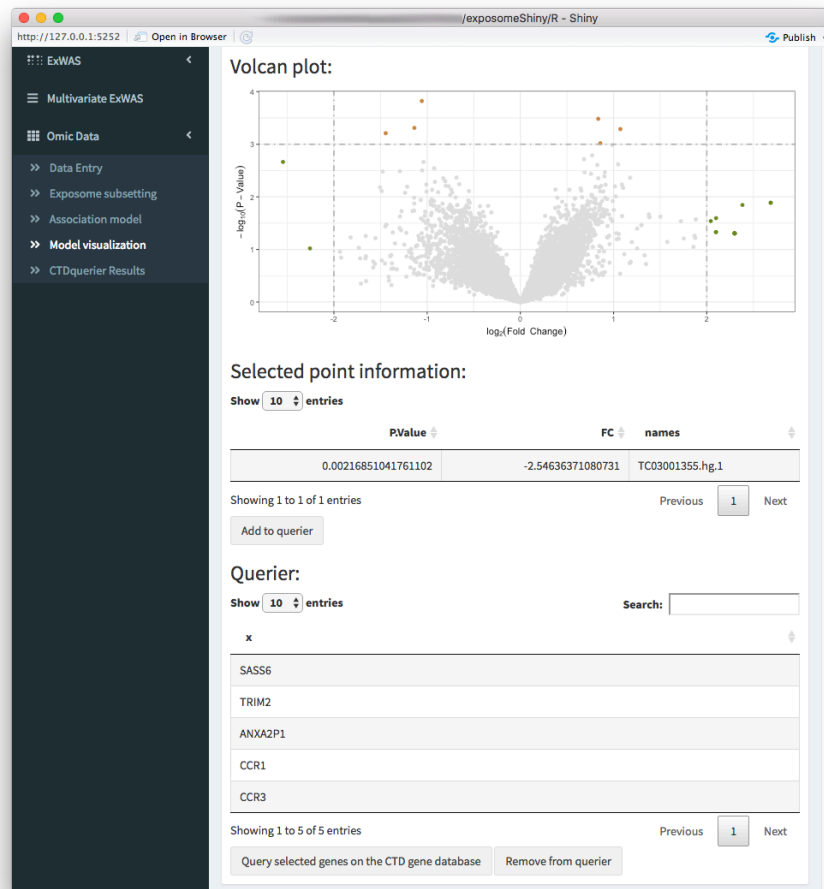

To visualize the results of the query, go to the “CTDquerier results” subtab. There are six tabs showing different results interpretations. First there’s the “Lost & found” tab which a plot to see the amount of genes found on the CTD database and the ones that were not found them, ther’s also two lists stating the names of them.

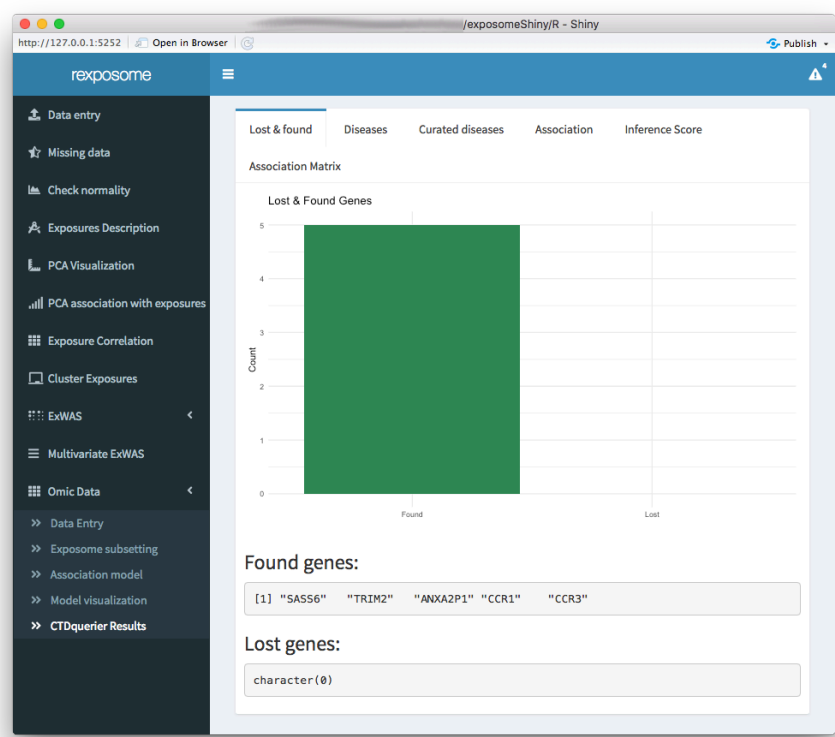

The diseases tab shows a table of all the associated diseases found on the CTD database.

	Disease.Name	Disease.ID	Direct.Evidence	Inference.Score	Reference.Count	GeneSymbol	GeneID
1	MICROCEPHALY 14, PRIMARY, AUTOSOMAL RECESSIVE	OMIM:616402	marker/mechanism		1	SASS6	163786
2	Weight Loss	MESH:D015431		91.69	80	SASS6	163786
3	Inflammation	MESH:D007249		68.53	85	SASS6	163786
4	Chemical and Drug Induced Liver Injury	MESH:D056486		68.35	568	SASS6	163786
5	Poisoning	MESH:D011041		66.54	29	SASS6	163786
6	Necrosis	MESH:D009336		65.72	201	SASS6	163786
7	Fibrosis	MESH:D005355		58.63	24	SASS6	163786
8	Prenatal Exposure Delayed Effects	MESH:D011297		54.78	198	SASS6	163786
9	Abnormalities, Drug-Induced	MESH:D000014		50.33	60	SASS6	163786
10	Lung Neoplasms	MESH:D008175		47.55	93	SASS6	163786

Showing 1 to 10 of 6,350 entries


Previous 1 2 3 4 5 ... 635 Next

The curated diseases tab shows the table of associated diseases but only shows the ones with direct evidence.


The screenshot shows the 'rexposome' web application interface. The left sidebar contains navigation options: Data entry, Missing data, Check normality, Exposures Description, PCA Visualization, PCA association with exposures, Exposure Correlation, Cluster Exposures, ExWAS, Multivariate ExWAS, and Omic Data. The main content area is titled 'Association' and displays a table of disease-gene associations. The table has columns for Disease.Name, Disease.ID, Inference.Score, Reference.Count, GeneSymbol, and GeneID. The table lists 10 entries, with the first entry being MICROCEPHALY 14, PRIMARY, AUTOSOMAL RECESSIVE (OMIM:616402) associated with the SASS6 gene (GeneID: 163786). The second entry is Osteoarthritis (MESH:D010003) associated with the TRIM2 gene (GeneID: 23321). The third entry is CHARCOT-MARIE-TOOTH DISEASE, AXONAL, TYPE 2R (OMIM:615490) associated with the TRIM2 gene (GeneID: 23321). The fourth entry is Liver Diseases (MESH:D008107) associated with the CCR1 gene (GeneID: 1230). The fifth entry is Carcinoma, Hepatocellular (MESH:D006528) associated with the CCR1 gene (GeneID: 1230). The sixth entry is Pneumonia (MESH:D011014) associated with the CCR1 gene (GeneID: 1230). The seventh entry is Hypersensitivity (MESH:D006967) associated with the CCR1 gene (GeneID: 1230). The eighth entry is Dermatitis, Contact (MESH:D003877) associated with the CCR1 gene (GeneID: 1230). The ninth entry is Behcet Syndrome (MESH:D001528) associated with the CCR1 gene (GeneID: 1230). The tenth entry is Respiratory Hypersensitivity (MESH:D012130) associated with the CCR3 gene (GeneID: 1232). The interface also includes a search bar, a 'Show 10 entries' dropdown, and pagination controls at the bottom.

	Disease.Name	Disease.ID	Inference.Score	Reference.Count	GeneSymbol	GeneID
1	MICROCEPHALY 14, PRIMARY, AUTOSOMAL RECESSIVE	OMIM:616402		1	SASS6	163786
2	Osteoarthritis	MESH:D010003	4.37	5	TRIM2	23321
3	CHARCOT-MARIE-TOOTH DISEASE, AXONAL, TYPE 2R	OMIM:615490		1	TRIM2	23321
4	Liver Diseases	MESH:D008107	84.12	77	CCR1	1230
5	Carcinoma, Hepatocellular	MESH:D006528	63.91	229	CCR1	1230
6	Pneumonia	MESH:D011014	56.49	52	CCR1	1230
7	Hypersensitivity	MESH:D006967	27.87	14	CCR1	1230
8	Dermatitis, Contact	MESH:D003877	17.96	14	CCR1	1230
9	Behcet Syndrome	MESH:D001528	2.43	4	CCR1	1230
10	Respiratory Hypersensitivity	MESH:D012130	32.38	11	CCR3	1232


The association tab shows information about all the direct evidence associated diseases. Select the disease of interest to see the score and reference count of it.

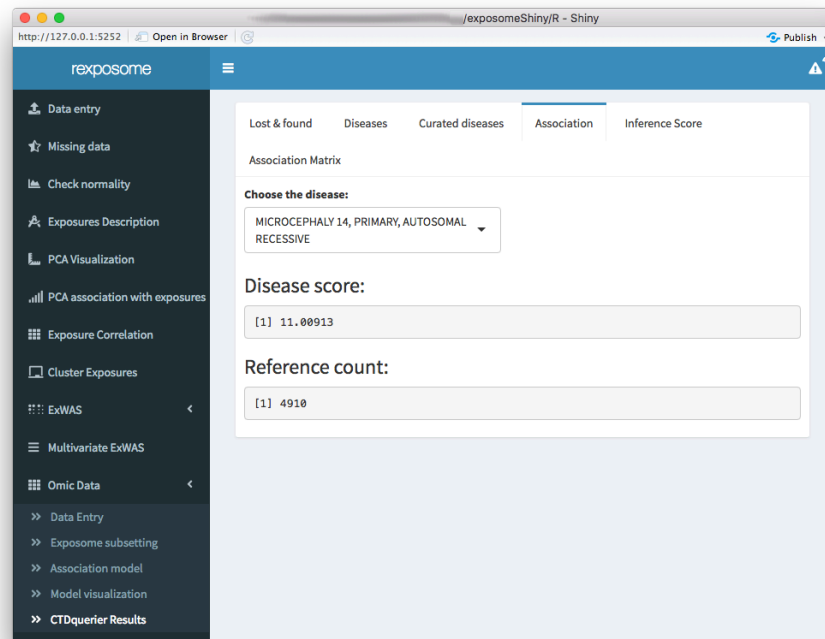

The inference score tab shows the inference score for each gene for a selected disease, the filter parameters puts out the genes with an inference score lower than the selected filter.

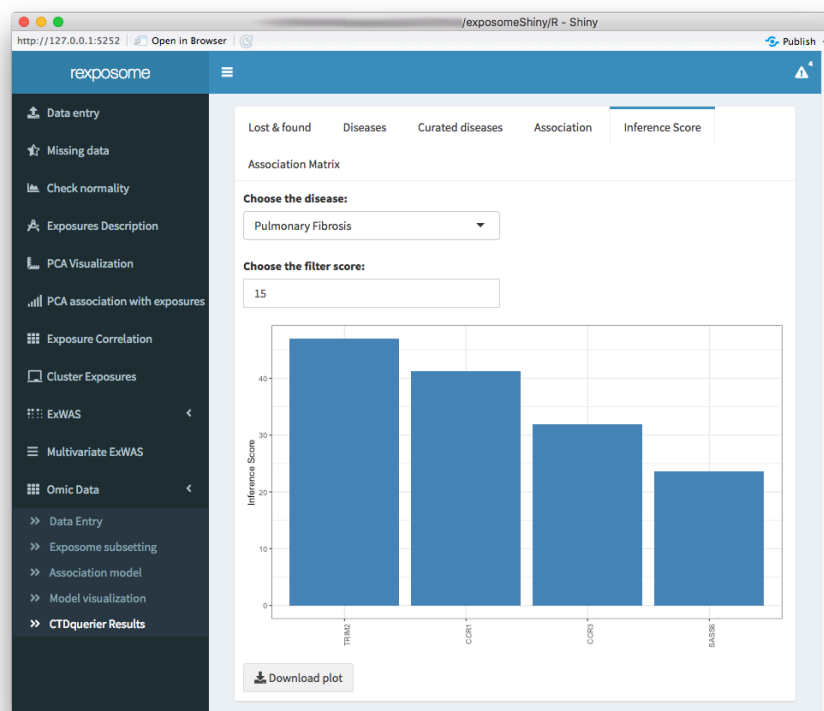


The association matrix tab shows a matrix of genes vs. chemicals with a heatmap representing the existing papers (references) providing evidence about the association between chemicals and genes.

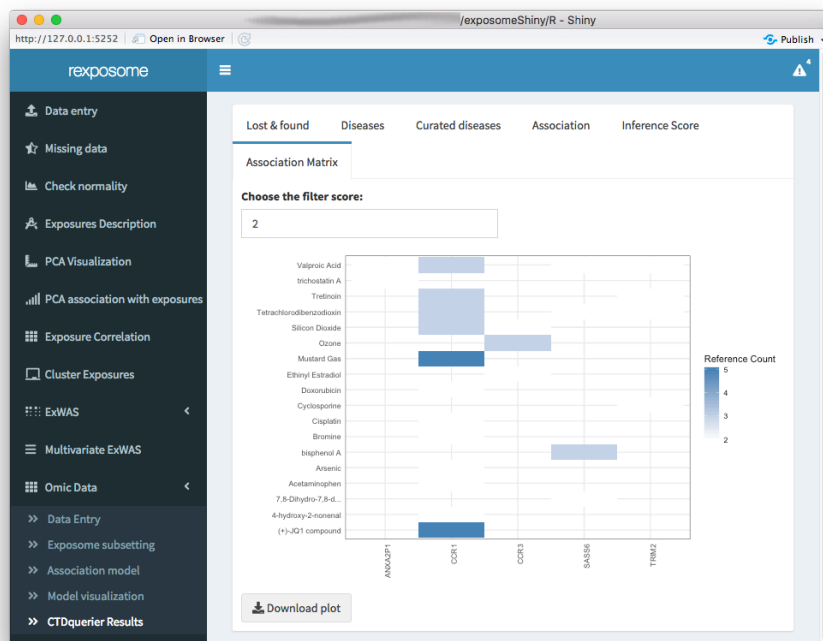

Chapter 6


General application functionalities

This whole Shiny application serves the purpose to perform a various number of exposome and omic analysis with an input set of data. To perform this analysis some operations can performed on the inputted dataset prior to the analysis and so in order to follow track of what exactly has been done or what is loaded on the current session, there's implemented some sort of state tracker inside the Shiny application. In order to access it, press the icon on the top right of the application

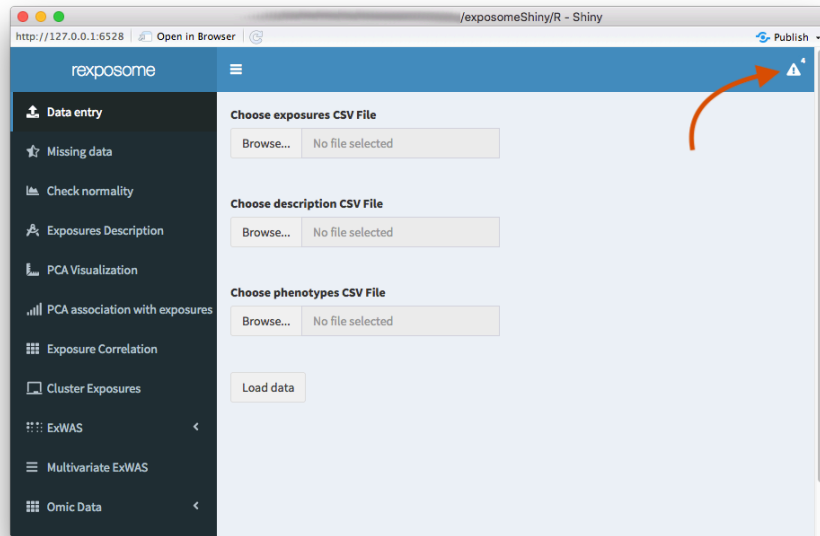

When clicking it a dropdown menu appears, inside there are seven different notifications:

- Exposome dataset: Turns to 100% when an exposome dataset is loaded in the environment. Here's a graphical example.

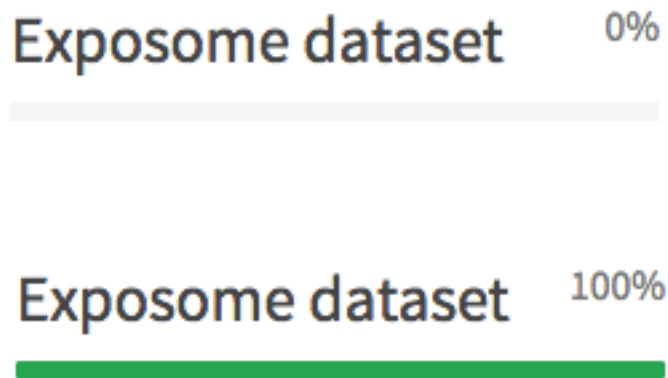

- LOD imputed: Turns to 100% if the exposome dataset is LOD imputed.
- Missing imputed: Turns to 100% if the exposome dataset has the missings imputed.
- Normality corrected: Turns to 100% if the exposome dataset is normality

corrected.

- Omics dataset: Turns to 100% when an omics dataset is loaded in the environment.
- Subset: Turns to 100% (and displays the subset family(ies)) when the exposome dataset is subseted.
- Model: Turns to 100% (and displays the association variable(s)) when a model is performed for an omics association analysis. Here's an example of a subset and model information.

Subset: Metals 100%



Model: cbmi 100%

